

PROPOR SG Filter Cartridges

- liquid filters
- polyethersulphone



Note: PROPOR and DEMICAP are registered trademarks of Parker domnick hunter



PROPOR SG sterilizing grade filters feature a patented, microbially retentive polyethersulphone membrane for fast, reliable and cost-effective sterile filtration of pharmaceutical fluids.


The asymmetric pore structure and high voids volume of the PROPOR SG membrane allow high throughputs and exceptionally high flow rates compared with competitive PES and alternative membranes. Low protein and preservative binding properties minimize product loss due to adsorption.

PROPOR SG filters are optimized for pharmaceutical processing. They have low extractable levels and broad chemical compatibility across the full pH range including organic solvents.


Features and Benefits

- Up to 3.5 times higher flow rates than competitive sterilizing grade filters
- Fully validated and integrity testable membrane for assurance of sterility
- Low binding for minimal product loss
- MURUS and DEMICAP's can be gamma-irradiated and autoclaved


Performance Characteristics

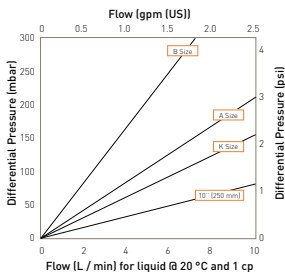

Cartridge flow rates
0.2 µm Cartridge

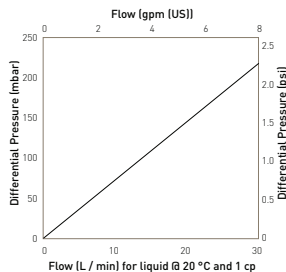

MURUS flow rates (10" Size [250 mm])
0.2 µm Capsule

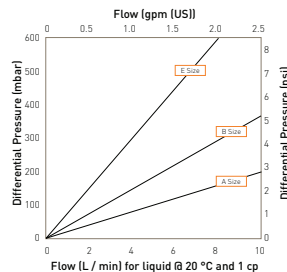

DEMICAP flow rates
0.2 µm Capsule

PROPOR SG Filter Cartridges

Specifications

Materials of Construction

- Filtration Membrane: Polyethersulphone
- Upstream Support: Polyester
- Downstream Support: Polyester

Filter Cartridges

- Inner Support Core: Polypropylene
- Outer Protection Cage: Polypropylene
- End Caps: Nylon
- End Caps Insert: 316L Stainless Steel

MURUS Disposable Filter Capsules

- Core: Polypropylene
- Sleeve: Polypropylene
- End Caps Insert: 316L Stainless Steel
- Standard o-rings/gaskets: Silicone
- Capsule Body: Polypropylene
- Capsules Vent Seals: Silicone

DEMICAP Filter Capsules

- Core: Polypropylene
- Sleeve: Polypropylene
- End Caps: Nylon
- Capsule Body: Nylon
- Capsules Vent Seals: Silicone
- Filling Bell: Polycarbonate

Syringe Filters

- Body: Polypropylene

Recommended Operating Conditions

Filter Cartridges

Up to 70 °C (158 °F) continuous operating temperature and higher short-term temperatures during CIP to the following limits:

Temperature °C	Temperature °F	Max. Forward dP (bar)	Max. Forward dP (psi)
20	68	5.0	72.5
40	104	4.0	58.0
60	140	3.0	43.5
80	176	2.0	29.0
90	194	1.7	24.6

MURUS Disposable Filter Capsules

Up to 25 °C (77 °F) @ 5.5 barg (79.7 psig)
Up to 60 °C (140 °F) @ 2.8 barg (40.6 psig)

Parker Hannifin certifies that this product complies with the European Council Pressure Equipment Directive (PED) 97/23/EC Article 3, Paragraph 3 - Sound Engineering Practice (SEP). This product is intended for use with Group 1 & 2 Dangerous and Harmless Liquids and Group 2 Harmless Gases at the operating conditions stated in this document. In compliance with PED Article 3, Paragraph 3, SEP, this product does not bear the CE mark.

DEMICAP Filter Capsules

Up to 40 °C (104 °F) at line pressures up to 5.0 barg (72 psig).

Effective Filtration Area (EFA)

10" (250 mm):	0.55 m ²	[5.92 ft ²]
K Size:	0.26 m ²	[2.79 ft ²]
A Size:	0.20 m ²	[2.15 ft ²]
B Size:	0.10 m ²	[1.07 ft ²]
E Size:	0.05 m ²	[0.53 ft ²]
Syringe ø50 mm:	14.50 cm ²	[2.25 in ²]

Sterilization

	Autoclave Cycles		Steam-in-Place Cycles (30 min.)	
	Cycles	Temp	Cycles	Temp
Cartridges	10	130 °C (266 °F)	30	130 °C (266 °F)
MURUS	5	130 °C (266 °F)	-	-
DEMICAP	10	130 °C (266 °F)	-	-
Syringe	1	130 °C (266 °F)	-	-

PROPOR SG filter cartridges can be sanitized with hot water at up to 90 °C (194 °F) and are compatible with a wide range of chemicals.

For detailed operational procedures and advice on cleaning and sterilization, please contact the Technical Support Group through your usual Parker domnick hunter contact.

Food and Biological Safety

Materials conform to the relevant requirements of 21CFR Part 177 and current USP Plastics Class VI - 121 °C and ISO10993 equivalents.

Quality Standards

Pharmaceutical grade products are manufactured in accordance with cGMP, 100% flushed with pharmaceutical grade purified water and integrity tested prior to despatch. A sample of each lot is tested to demonstrate conformity to validated claims.

Gamma-Irradiation

PROPOR SG MURUS & DEMICAP disposable filters can be gamma-irradiated up to a maximum dosage of 40 kGy.

Performance Characteristics

TOC / Conductivity
The filtrate quality from a 10" (250 mm) PROPOR SG conforms to the requirements of current USP <643> (TOC) and USP <645> (conductivity) within the first 200 ml flush of purified water.

Endotoxins
Aqueous extracts from the 10" (250 mm) PROPOR SG contain < 0.25 EU / ml when tested in accordance with the Limulus Amoebocyte Lysate test.

Non-Volatile Extractables (NVE)
Total NVEs extracted in the first 5 litre flush of purified water for a 10" (250 mm) cartridge are <10 mg.

Total NVEs extracted in the first 5 litre flush of purified water for an A size 7.9" (200 mm) DEMICAP capsule are <5 mg.

Pharmaceutical Validation
A full validation guide is available upon request from Laboratory Services Group (LSG).


Oxidizable Substances
PROPOR SG filter cartridges meet current USP and EP quality standards for sterile purified water for oxidizable substances following a <1 litre water flush.

Integrity Test Data
All filters are integrity testable to the following limits when wet with water and using air as the test gas.


Micron Rating	0.1	0.2	0.45
Filter Cartridges / MURUS / DEMICAP / Syringe Filters			
Min. Bubble Point* (barg)	2.36	3.38	2.48
(psig)	34.2	49.0	36.0
Filter Cartridges / MURUS / DEMICAP / Syringe Filters			
Diffusional Flow (barg)	4.8	2.8	1.7
Test Pressure (psig)	69.6	40.6	24.9
Filter Cartridges / MURUS / DEMICAP / Syringe Filters			
Max. Diffusional Flow (10 ⁻¹) (ml / min)	27.0	16.0	16.0
(K)	12.6	7.5	7.5
(A)	10.1	6.0	6.0
(B)	4.9	2.9	2.9
(E)	2.1	1.2	1.2

*Bubble point for 0.1 µm product is in 60/40 v/v IPA/Water.

Retention Characteristics
PROPOR SG filter cartridges are validated by bacterial challenge testing with *Brevundimonas diminuta* to current ASTM F838-05 methodology [10⁷ organisms / cm² EFA minimum] with typical in-house challenge levels being 10¹¹ organisms per 10" (250 mm) filter cartridge.

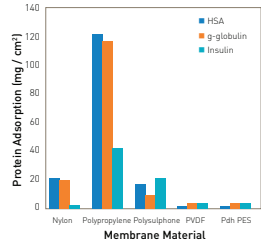

Protein binding on membrane materials

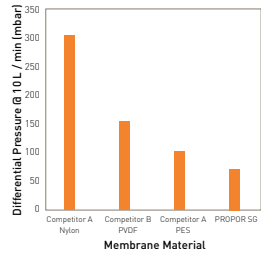

Differential pressure comparison of 10" (250 mm) sterilising grade filters

Ordering Information

Cartridges

ZCSG

Code Length (Nominal)	Code Micron	Code Endcap (10")	Code Variant	Code O-rings	Code Componentry
B* 2.5" (65 mm)	010 0.10 µm	C BF / 226 Bayonet	P Pharmaceutical	E EPDM	X All polypropylene componentry
A* 5" (125 mm)	020 0.20 µm	D Fin / 222		S Silicone	
K 5" (125 mm)	045 0.45 µm	E Flat Top / 222		V Viton	
1 10" (250 mm)		G Recess / 222			
2 20" (500 mm)		R BF / 222 Bayonet			
3 30" (750 mm)					
4 40" (1000 mm)					

* Supplied in packs of 3.

Code Endcap (Demi)
SK Retrofit
T TRUESEAL
Y Demi Stub
Z Demi A & B Std

MURUS Capsules

ZLSG

Code Length (Nominal)	Code Micron	Code Inlet Connection	Code Outlet Connection	Code Variant	Code Grade	Code Design	Code O-rings
K 5" (125 mm)	010 0.10 µm	A 1/2" Tri-Clamp	A 1/2" Tri-Clamp	P Pharmaceutical	N Non-sterile	L In-Line	E EPDM
1 10" (250 mm)	020 0.20 µm	B 1 1/2" Tri-Clamp	B 1 1/2" Tri-Clamp		S Pre-sterilized	T T-Port	V Silicone
2 20" (500 mm)	045 0.45 µm	D 1" Hosebarb	D 1" Hosebarb		γ (>25 kGy)	X* All polypropylene componentry	
3 30" (750 mm)		T 1" Tri-Clamp	T 1" Tri-Clamp				

*Only available with a 1" Tri-Clamp

*Silicone o-ring supplied as standard without flaring to specify the 'S' code.

DEMICAP Capsules

ZESG

Code Length (Nominal)	Code Micron	Code Inlet Connection	Code Outlet Connection	Code Variant	Code Grade	Code Pack N°	Code Accessory
E 4.4" (113 mm)	010 0.10 µm	T 1" Tri-Clamp	T 1" Tri-Clamp	P Pharmaceutical	N Non-sterile	3 Pack of 3	FB Filling Belt
B 5.5" (140 mm)	020 0.20 µm	N 1/2" NPT Male	N 1/2" NPT Male		S Pre-sterilized	X* All polypropylene componentry	
A 7.9" (200 mm)	045 0.45 µm	H 1/2" Hosebarb	H 1/2" Hosebarb		γ (>25 kGy)		
		G Stepped Hosebarb	G Stepped Hosebarb				
		M 1/2" NPT Male	M 1/2" NPT Male				
		Q Walther QC	Q Walther QC				
		R Grommel / QC	R Grommel / QC				

*Note: Supplied in Packs of 3

® & H connections only

Syringe Filters

ZSSG

Code Diameter	Code Micron	Code Inlet Connection	Code Outlet Connection	Code Variant	Code Grade	Code Options	Code Pack N°
050 50 mm	010 0.10 µm	F Female Luer Lock	F Female Luer Lock	P Pharmaceutical	N Non-sterile	S Standard	025 25 per box
	020 0.20 µm	G Stepped Hosebarb	G Stepped Hosebarb				
	045 0.45 µm						

Parker danmark herfor has a continuous policy of product development and although the Company reserves the right to change specifications, it attempts to keep customers informed of any alterations. This publication is for general information only and customers are requested to contact our Process Filtration Sales Department for detailed information and advice on a products suitability for specific applications. All products are sold subject to the company's Standard conditions of sale.